


# Advantages of the Super Grip-it

Designed for use in top-quality  
surgical instruments

Our process uses Creep Feed technology to diamond grind individual teeth into the hard, wear resistant tungsten carbide to produce precise, super sharp teeth with a well-defined point and inserts with superb gripping properties.



Diamond Grinding Wheel.

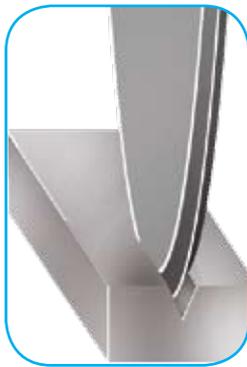


Diamond Grinding Wheel forms a "V" shape into the face of the carbide inserts.



The inserts are ground to create sharp grooves.



Diamond Grinding Wheel forms a "V" through the inserts to create the teeth.

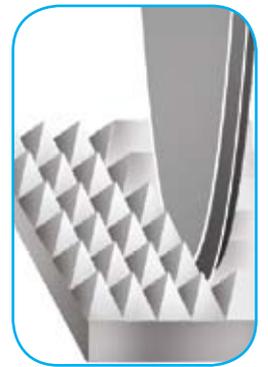

The finished result is super sharp, precise teeth across the inserts.

## Contact

We can be contacted in a number of ways:

By post:  
Serrations Ltd, 5 Cobnar Wood Close,  
Chesterfield Trading Estate,  
Chesterfield, United Kingdom S41 9RQ.

By telephone: +44 (0) 1246 456595

By fax: +44 (0) 1246 453377

By email:  
enquiries@serrations.com